


# 20075

Sterilisation packaging

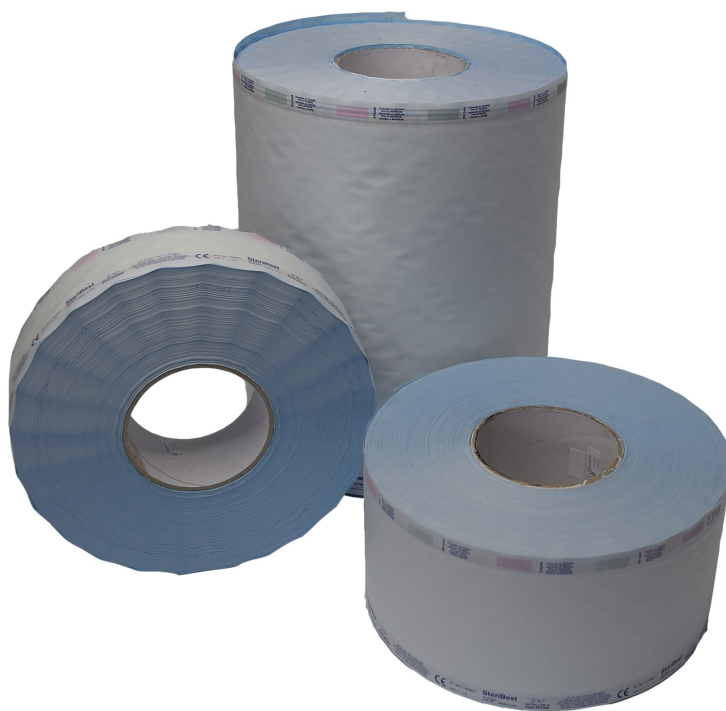

High tear and puncture resistance

Secure seal seams

Reliable sterile barrier

MADE IN  
EUROPE

## 20075

**Material:** Polypropylene / Paper, blue

**Miscellaneous:** Width of roll 75 mm

**Packaging unit:** 8 piece(s)

Sterilisation foil with steam and EO gas indicator,, blue (colour code: 2200), polypropylene /paper, high tear and puncture resistance, roll length 200 m, roll width 75 mm, suitable for F0 sterilisation, secure seal seams, reliable sterile barrier, sealing temperature: 160°C / 180°C, time between the machine's sealing bars and the product: 0.9 seconds / 1.5 seconds, pressure: 1–3 bar, recommended temperature: 160°C and 180°C, Steam sterilisation: 121°C = 15 minutes (minimum) / 134°C = 5 minutes (minimum), EtO gas: 50–60°C = 3–6 hours



### Industries

-  Beauty Edition
-  Black Edition
-  Dental Edition
-  Medical Edition

### Advantages

- + High tear and puncture resistance
- + Secure seal seams
- + Reliable sterile barrier



## Additional information

### Logistics data

| Size     | GTIN (SPU)    | GTIN (PU)     | Length (m) | Width (m) | Height (m) | Weight (kg) | Quantity PU / pallet |
|----------|---------------|---------------|------------|-----------|------------|-------------|----------------------|
| One size | 4059904060906 | 4059904029712 | 0,43 m     | 0,33 m    | 0,35 m     | 13,79 kg    | 30                   |

Customs tariff number: 48025525000

Country of origin: Italy