


20075

Sterilisation packaging

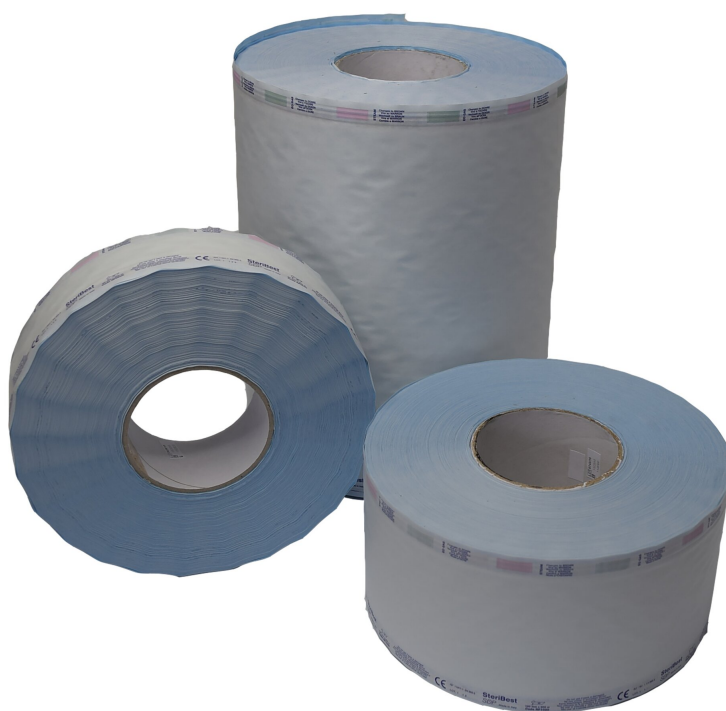

High tear and puncture resistance

Secure seal seams

Reliable sterile barrier

MADE IN
EUROPE

20075

Material: Polypropylene / Paper, blue

Miscellaneous: Width of roll 75 mm

Packaging unit: 8 piece(s)

Sterilisation foil with steam and EO gas indicator,, blue (colour code: 2200), polypropylene /paper, high tear and puncture resistance, roll length 200 m, roll width 75 mm, suitable for F0 sterilisation, secure seal seams, reliable sterile barrier, sealing temperature: 160°C / 180°C, time between the machine's sealing bars and the product: 0.9 seconds / 1.5 seconds, pressure: 1–3 bar, recommended temperature: 160°C and 180°C, Steam sterilisation: 121°C = 15 minutes (minimum) / 134°C = 5 minutes (minimum), EtO gas: 50–60°C = 3–6 hours



Industries

-  Beauty Edition
-  Black Edition
-  Dental Edition
-  Medical Edition

Advantages

- + High tear and puncture resistance
- + Secure seal seams
- + Reliable sterile barrier



Additional information

Logistics data

Size	GTIN (SPU)	GTIN (PU)	Length (m)	Width (m)	Height (m)	Weight (kg)	Quantity PU / pallet
One size	4059904060906	4059904029712	0,43 m	0,33 m	0,35 m	13,79 kg	30

Customs tariff number: 48025525000

Country of origin: Italy